

Product datasheet

info@arigobio.com

ARG56180 anti-Nucleoli Marker antibody [NM95]

Package: 50 μg Store at: -20°C

Summary

Product Description Mouse Monoclonal antibody [NM95] recognizes Nucleoli Marker

Tested Reactivity Hu

Tested Application FACS, ICC/IF, IHC-P

Host Mouse

Clonality Monoclonal

Clone NM95

Isotype IgG1, kappa

Target Name Nucleoli Marker

Species Human

Immunogen Nuclei of myeloid leukemia biopsy cells.

Conjugation Un-conjugated

Application Instructions

Application table	Application	Dilution
	FACS	0.5 - 1 μg/10^6 cells
	ICC/IF	0.5 - 1 μg/ml
	IHC-P	0.5 - 1 μg/ml
Application Note	IHC-P: Antigen Retrieval: Boiling tissue sections in 10mM citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 min. * The dilutions indicate recommended starting dilutions and the optimal dilutions or concentrations should be determined by the scientist.	

Properties

Form Liquid

Purification Purification with Protein G.

Buffer PBS (pH 7.4), 0.05% Sodium azide and 0.1 mg/ml BSA.

Preservative 0.05% Sodium azide

Stabilizer 0.1 mg/ml BSA


Concentration 0.2 mg/ml

Storage instruction For continuous use, store undiluted antibody at 2-8°C for up to a week. For long-term storage, aliquot

and store at -20°C or below. Storage in frost free freezers is not recommended. Avoid repeated freeze/thaw cycles. Suggest spin the vial prior to opening. The antibody solution should be gently mixed


before use.

Images

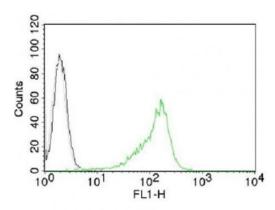


ARG56180 anti-Nucleoli Marker antibody [NM95] FACS image

Flow Cytometry: 293 cells stained with ARG56180 anti-Nucleoli Marker antibody [NM95] (green).

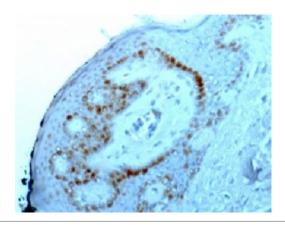


ARG56180 anti-Nucleoli Marker antibody [NM95] IHC-P image

Immunohistochemistry: Formalin-fixed and paraffin-embedded Human skin stained with ARG56180 anti-Nucleoli Marker antibody [NM95].

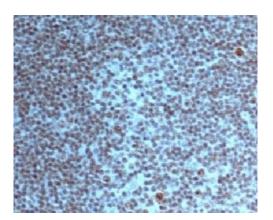

ARG56180 anti-Nucleoli Marker antibody [NM95] IHC-P image

Immunohistochemistry: Formalin-fixed and paraffin-embedded Human tonsil tissue stained with ARG56180 anti-Nucleoli Marker antibody [NM95].

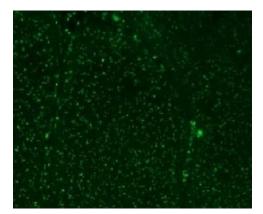

ARG56180 anti-Nucleoli Marker antibody [NM95] IHC-P image

Immunohistochemistry: Formalin-fixed and paraffin-embedded colon carcinoma stained with Alexa Fluor 488-conjugated ARG56180 anti-Nucleoli Marker antibody [NM95].