

Product datasheet

info@arigobio.com

ARG82184 Package: 100 tests
Monoamine Oxidase Inhibitor Screening Kit (Fluorometric) Store at: -20°C

Summary

Product Description ARG82184 Monoamine Oxidase Inhibitor Screening Kit (Fluorometric) is a screening kit for inhibitor

screening and evaluation of monoamine oxidase inhibitors. (Please note: the kit does not provide the

enzyme MAO-A and MAO-B.)

Tested Reactivity Hu, Ms, Rat, All

Tested Application FuncSt

Target Name Monoamine Oxidase Inhibitor

Conjugation Note Read at 530/585 nm (fluorometric)

Sample Volume 45 µl

Application Instructions

Application Note Please note that this kit does not include a microplate.

Assay Time 35 min

Properties


Form Liquid

Storage instruction Store the kit at -20°C. Do not expose test reagents to heat, sun or strong light during storage and usage.

Please refer to the product user manual for detail temperatures of the components.

Note For laboratory research only, not for drug, diagnostic or other use.

Images

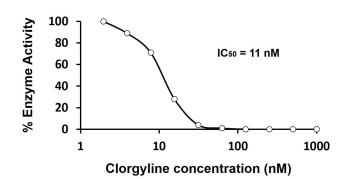

ARG82184 Monoamine Oxidase Inhibitor Screening Kit (Fluorometric) inhibitor titration curve image

Inhibitor titrations: Human MAO-A was incubated with various concentrations of Clorgyline. Each concentration of inhibitor contained 10 v/v% DMSO (final 0.5 v/v%). The IC50 for Clorgyline with 3 U/ml human MAO-A was determined to be 11 nM.